- Model name : MT-1000 / MT-3000
- Product name : Infrared skin thermometer
- Manufacturer : T&R BIOFAB
- Address : No. 545, Industry Convergence Head Office,  
Smart Hub, 237, Sangidaehak-ro, Siheung City, Gyeonggi-do
- Manufacture license no. 4627
- Manufacture certification no. : No. 16- 4397
- No of Mfg.(S/N) : Indicated on the product
- Date of Mfg. : Indicated on the product
- Packaging unit : 1 set

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
Product name	THERMOCARE Smart Thermometer
Date of purchase	
Place of purchase	
Model name	<input type="checkbox"/> MT-1000 <input type="checkbox"/> MT-3000
Warranty period	Within 1 year from the date of purchase

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
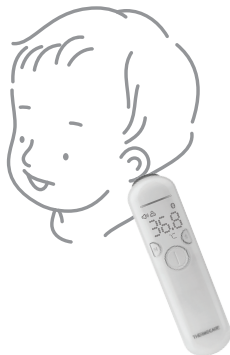
MEMO

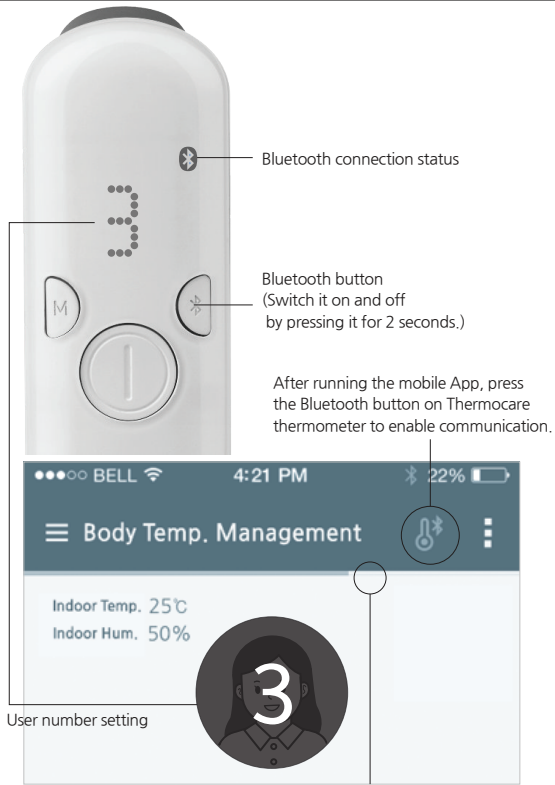The mobile app will discover the Thermocare thermometer registered on the device and automatically begin to fetch the measured body temperature information.

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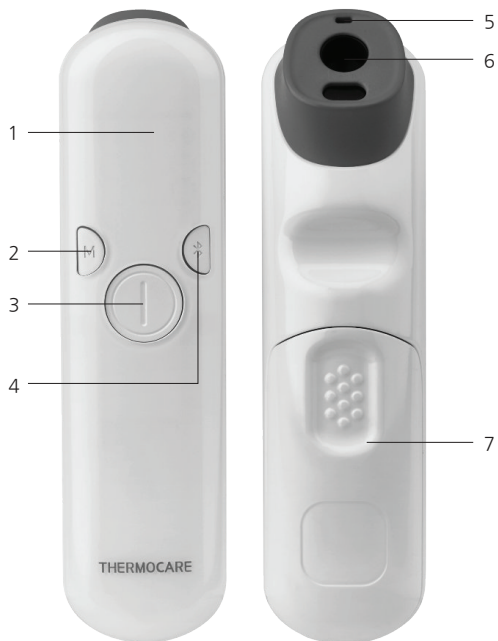
1. Install Thermocare App in the smart phone.

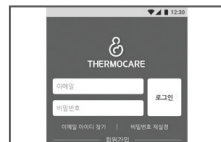2. After signing up, log-in.

3. Press 'Register device' at the top of the menu.

4. Turn on the Bluetooth function on the Thermocare thermometer. (If the Bluetooth button is pressed, the Bluetooth icon lights up).

5. Wait for a while. (Maximum 30 seconds)

- ※ Check if the Bluetooth in the smart phone has been turned on as well!
- ※ Only temperature that has been measured after Bluetooth connection has been achieved after battery has been inserted is sent to the mobile App.

THERMOCARE Smart Thermometer works with mobile app to manage body temperature data conveniently.

Android


Please search and install "THERMOCARE" in the Google Play Store

iPhone

Please search and install "THERMOCARE" in the Apple App Store

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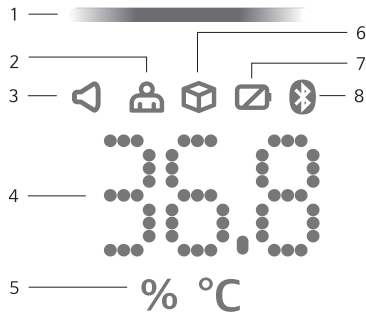How to keep and manage after use

- Keep out of water.
- Remove batteries when not in use for a long period of time.
- Stable state from vibration and impact.
- Keep out of the reach of children under the age of 12.
- Avoid fire and direct sunlight and store in a place free from dust and pollutants.
- If you find any signs of damage, contact Customer Center without attempting to repair yourself.
- It is recommended to store at normal room temperature, if stored outside room temperature, store the product at room temperature for minutes before use.
- The probe (temperature measurement part) of the thermometer is very sensitive, so keep it clean and clean in order to ensure accurate measurement.
- Keep it in a dry place.

Precautions when using

- Body temperature measured on the right forehead and on the left forehead may be different.
- Do not use this unit for any purpose other than its intended use.
- Measurement may be inaccurate when battery power is low. In this case, use after replacing the battery.
- Avoid touching the probe (temperature measurement part) of the thermocare thermometer or the battery part with your mouth by sucking or wet hands.
- Use after 30 minutes or more after the infant enters the room from outside or moves the thermometer to a different temperature environment.
- If you are sweating on the forehead, measure the area behind the earlobe rather than your forehead. After the hair around the ear is completely tilted, measure the probe (body temperature measuring part) against the back of the earlobe.
- Repeat three or more times in the following situations


- Infants within 90 days of age.
- If the infants is under 3 years of age and the measure temperature is out of the average temperature range.
- If you are new to the product or have no experience using similar products.

Precautions before use


- Learn how to use the product by reading the User Manual thoroughly.
- Use in surrounding temperature of 16°C ~ 40°C.  
(Use after it has been indoors for 30 minutes).
- Maintain an appropriate distance from direct fire that could interfere with the measurement.
- Check the connection status of the battery.
- Check if there is damage or abnormality in the device.
- Be aware of damage caused by impact.
- The body temperature may differ slightly (0.1~0.4°C) depending on the measurement area (mouth, ear, forehead).
- To save the battery, it is better to connect to Thermocare mobile App through Bluetooth after several measurements.

Operating condition and storage environment

- Operating temperature
- Storage temperature and humidity

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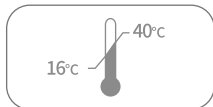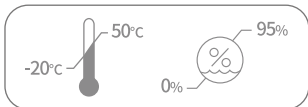
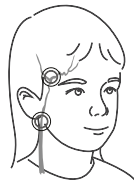
Object temperature measurement

- 1) Press the mode setting button to set the "Object, Sound On" state.
- 2) Press and release the measurement button located on the measurement object.
- 3) The measured temperature can be checked on the measurement display.
- 4) The object temperature measurement range 10°C to 70°C.

Measuring surrounding temperature/humidity

- 1) Press the Mode setting button for more than 3 seconds to measure the surrounding temperature/humidity.
  - 2) After the surrounding humidity is displayed for 3 seconds, the surrounding temperature is displayed.
- (※ Information on the surrounding temperature and humidity that have been measured must only be checked for reference.)
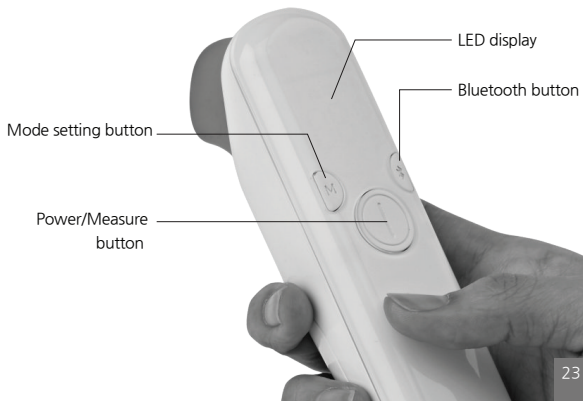
## Power

- 1) Insert 2 1.5V AA-type batteries into the product.
- 2) Press the Power/Measure button to turn on the power.

## Body Temperature measurement

- 1) Press the mode setting button to select one of "Human, Sound ON" / "Human, Sound OFF" / "Object, Sound ON".
- 2) Position the measurement site (temporal artery) and press the measurement button.
- 3) The measured body temperature can be checked on the measurement display and the high temperature/low temperature condition appears on the LED indicators.
- 4) When measuring body temperature, ambient temperature and humidity are measured together and can be confirmed by linking with mobile app.

- ① 36.0°C or less : orange light
- ② 36.0°C ~ 37.6°C : green light
- ③ 37.6°C ~ 38.5°C : orange light
- ④ exceeding 38.5°C : red light

- [Low temperature] Orange
- [Normal] Green
- [Slight fever] Orange
- [High temperature] Red

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

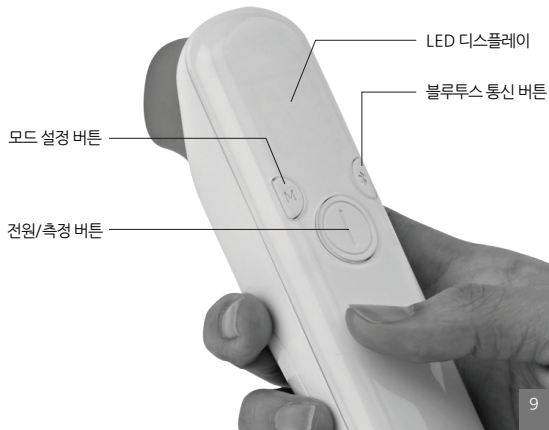
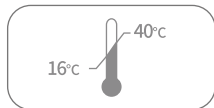
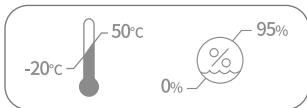
Purpose of use

Infrared skin thermometer is a device that measures body temperature by detecting the surface infrared emission of temporal arteries that is directly connected to the heart.

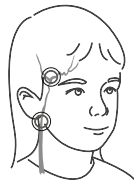Features

Name	Description
Rated voltage	DC 3V
Display method	LED display (measurement sound)
Measurement range	A. Body temperature: 22.0℃ ~ 42.2℃ B. Object temperature: 10.0℃ ~ 70.0℃ C. Humidity: 0%~100%
Accuracy	A. Body temperature: ± 0.3℃ (However, ± 0.2℃ in the range of 36.0℃ ~ 39.0℃) B. Object temperature: ± 2.0℃
Minimum scale unit	0.1 for temperature / 1 for humidity
Measurement time	Within 1 second
Memory capacity	1,000

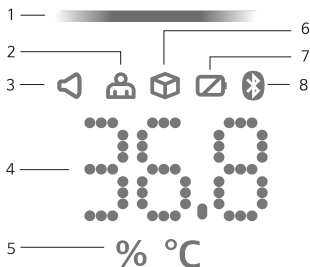No.	Name	Description
1	LED status display light	The color is displayed according to the temperature measurement range.
2	Temperature measurement mode	Displays the temperature measurement mode.
3	Measured sound indicator	Indicates sound ON/OFF.
4	Measuring temperature	The measured temperature is displayed. After the bluetooth connection, the selected measurement member number is displayed in the body temperature measurement mode.
5	Unit display	In case the surrounding humidity is indicated, % is displayed and as for body temperature, object temperature and surrounding temperature, °C is displayed.
6	Object measurement mode	Displays the object measurement mode.
7	Battery warning light	Below 10%, the warning light is displayed in red.
8	Bluetooth connection	Press the Bluetooth button until the LED on the front starts blinking blue.

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
No.	Name	Description
1	Measurement display	Measurement temperature, measurement mode, battery status
2	Measurement mode selection	Used when you select one of three modes. 1) Human, Measured sound ON 2) Human, Measured sound OFF 3) Object, Measured sound ON In case the button is pressed for 3 seconds, it displays the surrounding temperature/humidity information. (After displaying the surrounding humidity for 3 seconds, it displays the surrounding temperature)
3	Power/Measure button	It turns the power ON or measures the temperature.
4	Bluetooth button	Press and hold the button for 2 seconds to turn Bluetooth ON/OFF. (*After the Bluetooth connection, you can change the measurement member number by briefly pressing the button in the body temperature measurement mode.)
5	Temperature and humidity measuring sensor	Detects the temperature and humidity of the measurement environment.
6	Temperature measuring sensor	Detects the infrared rays emitted from the human body or object.
7	Battery	It contains 2 Alkaline AA batteries.

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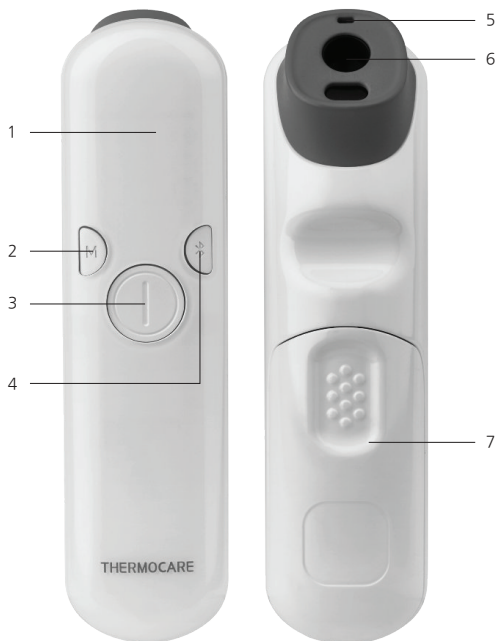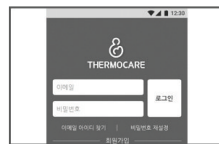
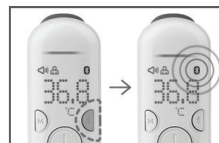
Table of contents

18-19	Name of each part
20	Screen description
21	Purpose of use and features
22-23	How to use
24	Precautions before use, operating conditions and storage environment
25	How to manage your storage and how to use it
26	Information on App download
27	Registering device
28	Bluetooth connection
29	Information on product warranty, MEMO

Normal times

1~2 cm above the eyebrows

It is based on the principle of measuring the temperature of blood flowing through the "temporal arteries" that passes by from the ear to the forehead.

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MEMO

When sweating

Area behind the earlobe to be measured

When there is sweat on the forehead, wipe the sweat before measuring or measuring the area behind the earlobe is recommended.

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# THERM&CARE

User Manual of Smart Thermometer, Thermocare